



# Step Ahead... IV Administration System

## SAFE...

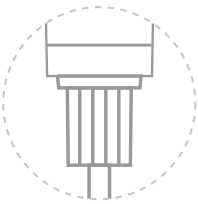


Eliminates chance of Air Embolism by **Air Stop**

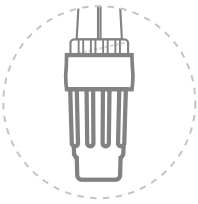

Ensures complete close system with **Prime Stop**



Neck Ring in Drip Chamber provides firm grip & increase paramedics convenience

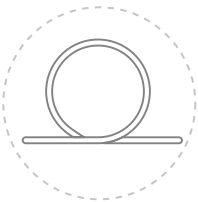

DEHP Free IV tube & chamber protects against phthalate toxicity

## STERILE...

- 0.2µ filter Air Vent maintains drip chamber sterility
- Hydrophilic **Air Stop** 15µ Filter provides efficient filtration
- **Prime Stop** prevents spillage and improves hygiene condition
- **Prime Stop** keeps the sterility intact till **Point of Administration**

# AeroVein™

## SUPERIOR...

Component	Advantage
Large Drip Chamber	<ul style="list-style-type: none"><li>■ Even drop formation and easy visual inspection</li><li>■ Neck ring provide firm grip while spiking the container</li><li>■ Lower soft chamber helps in effortless priming</li></ul>
Advance Roller Clamp	<ul style="list-style-type: none"><li>■ Provides accurate flow control</li></ul>
IV Tube	<ul style="list-style-type: none"><li>■ 180cm long DEHP free tube</li><li>■ Kink resistant to avoid tangling and maintain uninterrupted flow</li></ul>
Luer Lock	<ul style="list-style-type: none"><li>■ Safe and easy connection to cannula</li><li>■ Increases patient compliance</li><li>■ Eliminates risk of Phlebitis and Cellulitis</li></ul>
Spike Notch	<ul style="list-style-type: none"><li>■ Safe Disposal post usage</li></ul>
Tube Dock	<ul style="list-style-type: none"><li>■ To dock the tube during standby mode</li></ul>

# Step Ahead... IV Administration System

Feature	Benefits
Air Stop	<ul style="list-style-type: none"><li>Hydrophilic filter eliminates risk of air embolism</li></ul>
Prime Stop	<ul style="list-style-type: none"><li>Maintains close system</li><li>Improves hygiene standard</li><li>Sterility up to POA</li></ul>
Large Drip Chamber	<ul style="list-style-type: none"><li>Easy visual inspection</li><li>Effortless priming</li></ul>
Neck Ring	<ul style="list-style-type: none"><li>Firm grip while spiking</li></ul>
Advance Roller Clamp	<ul style="list-style-type: none"><li>Accurate flow control</li></ul>
Spike Notch	<ul style="list-style-type: none"><li>Safe Disposal</li></ul>
Tube Dock	<ul style="list-style-type: none"><li>To hold tubing</li></ul>
Lure Lock	<ul style="list-style-type: none"><li>Better patient compliance</li></ul>
0.2µ AirVent	<ul style="list-style-type: none"><li>Drip chamber Sterility</li></ul>



**AMANTA HEALTHCARE LTD.**

8<sup>th</sup> Floor, Shaligram Corporates,  
C. J. Marg, Ambli, Ahmedabad 380058, Gujarat, INDIA.  
Email : [info@amanta.co.in](mailto:info@amanta.co.in)  
Website: [www.amanta.co.in](http://www.amanta.co.in)

# AeroVein™

Vented I. V. Set with  
Air Stop & Prime Stop

